



SUPER•SAL™ UNIVERSAL SALIVA COLLECTION KIT, CATALOG # SSAL-601

INSTRUCTIONS FOR SALIVA COLLECTION USING THE SUPER•SAL™ UNIVERSAL SALIVA COLLECTION KIT


1. Do not eat, drink, smoke, or use oral hygiene products for at least 30 minutes before you start the collection process.
2. Prior to collection, review the detailed instructions for the Super•SAL™ Universal Saliva Collection Kit below, particularly the illustrations showing the appearance of the Sample Volume Adequacy Indicator (SVAI) changing from white to **fully red** after sample collection (Figure 4 below).
3. Open the package containing the Super•SAL™ Universal Saliva Collection Kit and place contents on a clean and dry surface. Follow the illustrated instructions below.

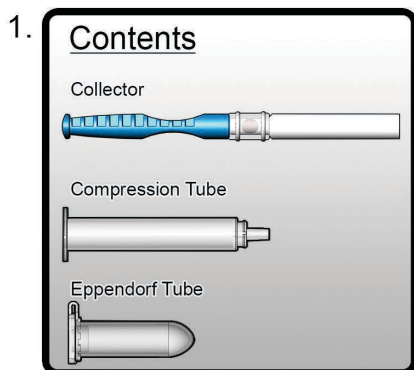


Place contents on a clean and dry surface. In preparation for saliva collection, **pool saliva** in the mouth. Video instructions are available on our website at: 4saliva.com/products/super-sal/

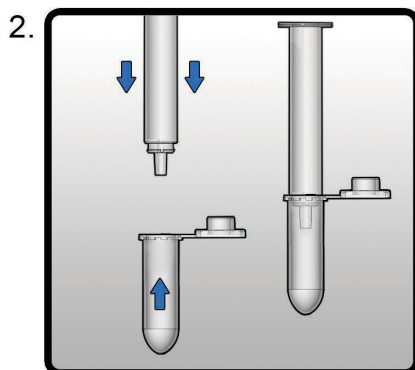


Attach the Eppendorf tube provided to the base of the Compression Tube firmly. Sample collection may now begin.

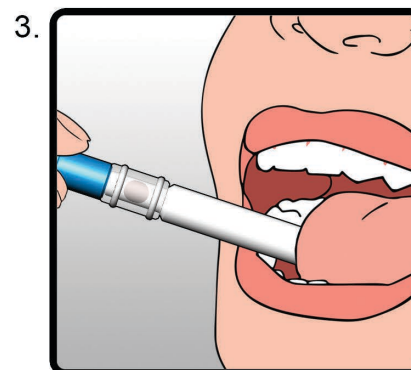


Place the tip of the pad of the Collector where saliva has pooled. The Collector may be removed from the mouth periodically to read the SVAI, but resume collection immediately afterward.

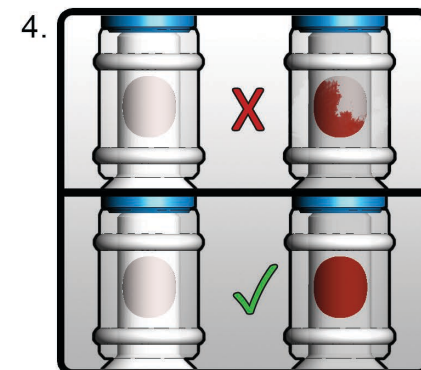


Collect until the appearance of the SVAI changes **completely to red**. The collection procedure should last approximately 1-5 minutes.

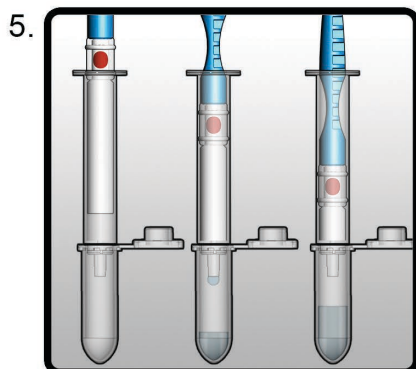

Place the white absorbent pad end into the Compression Tube holding the Collector in an upright and vertical position and firmly push the plunger downwards to transfer saliva from the absorbent pad into the Eppendorf Tube. Push and hold for 30 seconds.

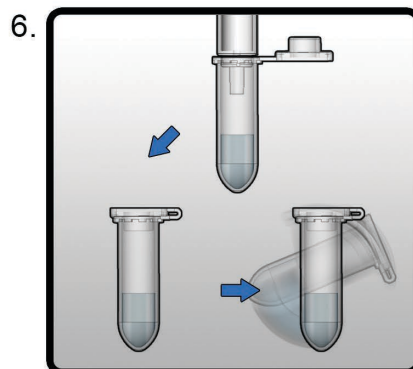

Gently remove the Eppendorf Tube from the end of the Compression Tube (it is appropriate to remove the Eppendorf Tube during expression if necessary due to foam blockage). Add preservative (not provided), if required, and close the lid on the tube tightly.

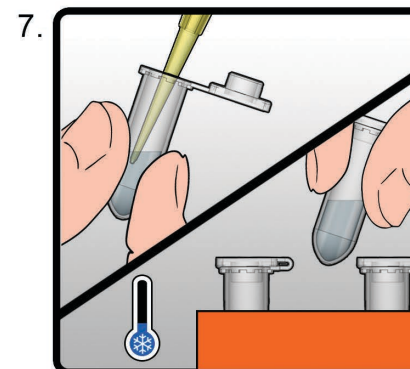

Process the collected sample immediately or store at refrigerated temperatures (depending upon the specific application) until required. For Technical Support, contact Oasis Diagnostics® at (360) 546-1563 or by e-mail to info@4saliva.com.